



# CAM SQUARED

M2 meter  
**The smart one**

Compact  
Alignment-Free  
Ultra short measurement cycle



# CAM SQUARED +

**A great choice for almost any lab or industrial application, the CAM SQUARED is Imagine Optic's innovative answer to the need for laser quality testing and  $M^2$  measurement.**

Finally an  $M^2$  meter as easy and quick to set up as a beam profiler.

## APPLICATIONS

Laser beam quality testing is of utmost importance in many laser-based applications where beam waist and beam divergence matter:

- + Manufacturing, machining, welding for fluence
- + Imaging, for resolution
- + Fiber optics, for coupling
- + Free space optical communications and laser radar systems (LIDAR) for better propagation through turbulent atmosphere.

CAM SQUARED performs multiple measurements :  $M^2$ , divergence, focus diameter, waist position, Rayleigh length, thermal effects.

## FEATURES

- + **ISO 11146 standard compliant.** The measurement of intensity combined with phase allows to generate 10 to  $\infty$  of intensity frames from which is calculated the  $M^2$  factor, such as described in the ISO 11146 standard.
- + **Self aligned.** CAM SQUARED requires no alignment, making setup quick and easy.
- + **Short measurement cycle.** CAM SQUARED requires no translation, making measurement cycle very short and the solution perfectly adapted to pulsed lasers and dynamic applications.
- + **Optics free.** As no mirrors nor lenses are necessary, there are no optics introducing aberrations which impair the beam quality. There are also no coatings limiting the range of use of the sensor.
- + **SM1 thread** on the front of the sensor for easy mounting of optical densities in order to adapt to the power of the laser to be tested.



# SPECIFICATIONS\*

## OPERATING SPECS

	M	L	XL / XXL (on request)	SWIR
Aperture dimensions	4.5 x 3.7 mm <sup>2</sup>	6.9 x 5.1 mm <sup>2</sup>	13.8 x 10.2 / 22 x 22 mm <sup>2</sup>	9.3 x 7.4 mm <sup>2</sup>
Recommended beam diameter min. @ 1/e <sup>2</sup> (min. @ 1/e <sup>3</sup> )	0.7 mm (0.8 mm)	0.7 mm (0.8 mm)	0.7 mm / 1.2 mm (0.8 mm / 1.4 mm)	1.6 mm (1.9 mm)
max. @ 1/e <sup>2</sup> (max. @ 1/e <sup>3</sup> )	3 mm (3.6 mm)	4.2 mm (5 mm)	8.2 mm / 17.8 mm (10.1 mm / 21.9 mm)	6 mm (7.3 mm)
Maximum acquisition frequency	125 Hz (USB 3.0) 30 Hz (GigE)	55 Hz (USB 3.0) 30 Hz (GigE)	30 Hz (USB 3.0) / 10 Hz (10GigE)	150 Hz (USB 3.0) 49 Hz (GigE)
Wavelength range	350 - 1100 nm			900 - 1700 nm
Minimum power	0.15 nW		0.15 nW / 0.7 nW	0.3 pW
External trigger	TTL signal			
Operating system	Windows 10 & 11			
Measurement cycle time	~ ms typically, depending on settings			
Travel range	not limited by translation stage			
Typical M <sup>2</sup> accuracy	5%			
Pulsed sources	full compatibility			
Damage thresholds	50 mW / cm <sup>2</sup> in CW mode 50 uJ / cm <sup>2</sup> in Pulsed mode			

## MISC

Dimensions (Height x Width x Length)	50 x 50 x 55 mm <sup>3</sup>	56 x 56 x 60 mm <sup>3</sup> / TBD	70 x 70 x 71 mm <sup>3</sup>
Weight for USB version	200 g	200 g / 800 g	250 g
Mounting configuration	horizontal or vertical		
Working temperature	15 - 30 °C		
Interface	USB 3.0 or optional GigE converter	USB 3.0 / 10GigE	USB 3.0 or option. GigE converter
Power consumption	3.1 W	3.6 W / 14 W	5 W

## OPTION

CAM SQUARED can be upgraded in option for wavefront sensing. In this case, in addition to the M<sup>2</sup> meter, you get access to a complete wavefront sensor with the following features (see HASO datasheets for more information):


Repeatability: <  $\lambda/200$  RMS

Absolute wavefront measurement accuracy: ~  $\lambda/100$  RMS

Wavefront error measurement provides detailed quantitative knowledge of the cause of aberrations and beam quality.

\*Subject to changes without further notice


## DIMENSIONS\*\* (mm)

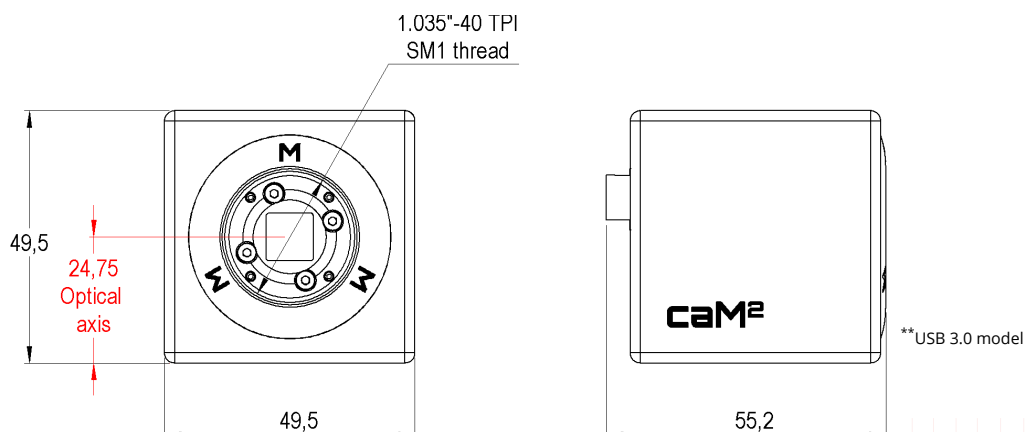

## SOFTWARE

### WAVESQUARED

- + Optimized display of laser quality metrics
- + Beam pointing adjustment and stability measurement
- + Optional phase measurement extension for wavefront diagnostic and analysis (alignment, collimation, optical aberrations analysis and more than 150 features)
- + Optional SDK in C/C++, LabVIEW and Python

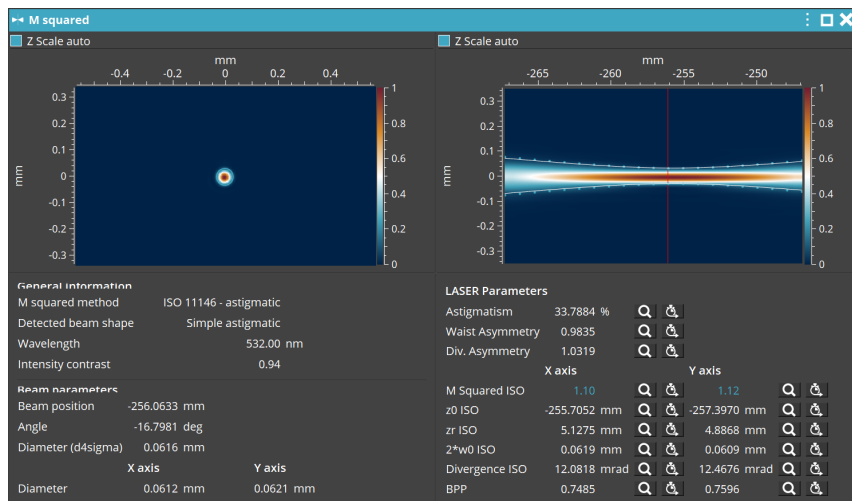

## ACCESSORIES

- + Several mounting options are available, including adapters for the most common mechanical stages and magnetically coupled top and bottom plates, allowing to mount, remove, and replace CAM SQUARED with a high repeatability.

## APPLICATION NOTES

- + M2 measurement with CAM SQUARED



## CONTACT US

**Imagine Optic Headquarters**  
18, rue Charles de Gaulle  
91400 ORSAY · France  
Phone +33 (0)1 64 86 15 60  
sales@imagine-optic.com  
www.imagine-optic.com

